

SC-100

COMPACT DENTAL STERILIZATION CENTER WATER PROCESSOR



The Steri-Center 100 is an ultra-compact water system perfect for your sterilization center. It simultaneously handles dental unit bottles, autoclaves, sterilizers and one instrument washer.

Features:

- Electric pump
 - **overcomes low tap water pressure**
 - consistent production at 4 gallons / hour
 - minimizes filter changes
- Produces 2 water types:
 - Mineral-free water for autoclave & sterilizer
 - Low-mineral water for dental unit bottles & instrument washer.
- Slim design for easy install in compact spaces
- Easy to maintain quick-change filters
- Includes
 - remote DI dispensing gun for direct fill of autoclave & sterilizer reservoirs
 - chrome faucet for dental unit bottle fill
 - 3.2 gallon storage tank

Advantages:

- **Minimizes operating disruption**
- **Eliminates angry dentists / product returns**
- Quick-change filters **eliminate service calls**
- **Filters last entire year** for most practices
- No mess QC filters keep your cabinet dry
- 1 year warranty


Benefits:

- Simultaneously handles autoclaves, sterilizers, dental unit bottles and an instrument washer.
- Eliminates bottled water hassles
 - Trips to purchase
 - Storage area
 - Lifting of heavy bottles
 - Staff injuries.
- Optional 6-gallon tank (DWS-ROSTC-10)
- Optional instrument washer kit (SC-1IWK)
 - 6-gallon tank
 - Miele & Hydrim hookup fittings

Specifications:

Tap water type	potable, unsoftened
Tap water temperature	40°F - 90°F
Tap water pressure	20 - 75 psi
Production	4 gph / 100 gpd
Height	16 in
Width	12 in
Depth	5 in
Water Ports	1/4"
Electrical (1 plug)	110VAC / 60Hz / 2A
Shipping weight	23 lbs
Annual service kit	SC-1PM
6-gallon tank	DWS-ROSTC-10
Instrument washer kit	SC-1IWK



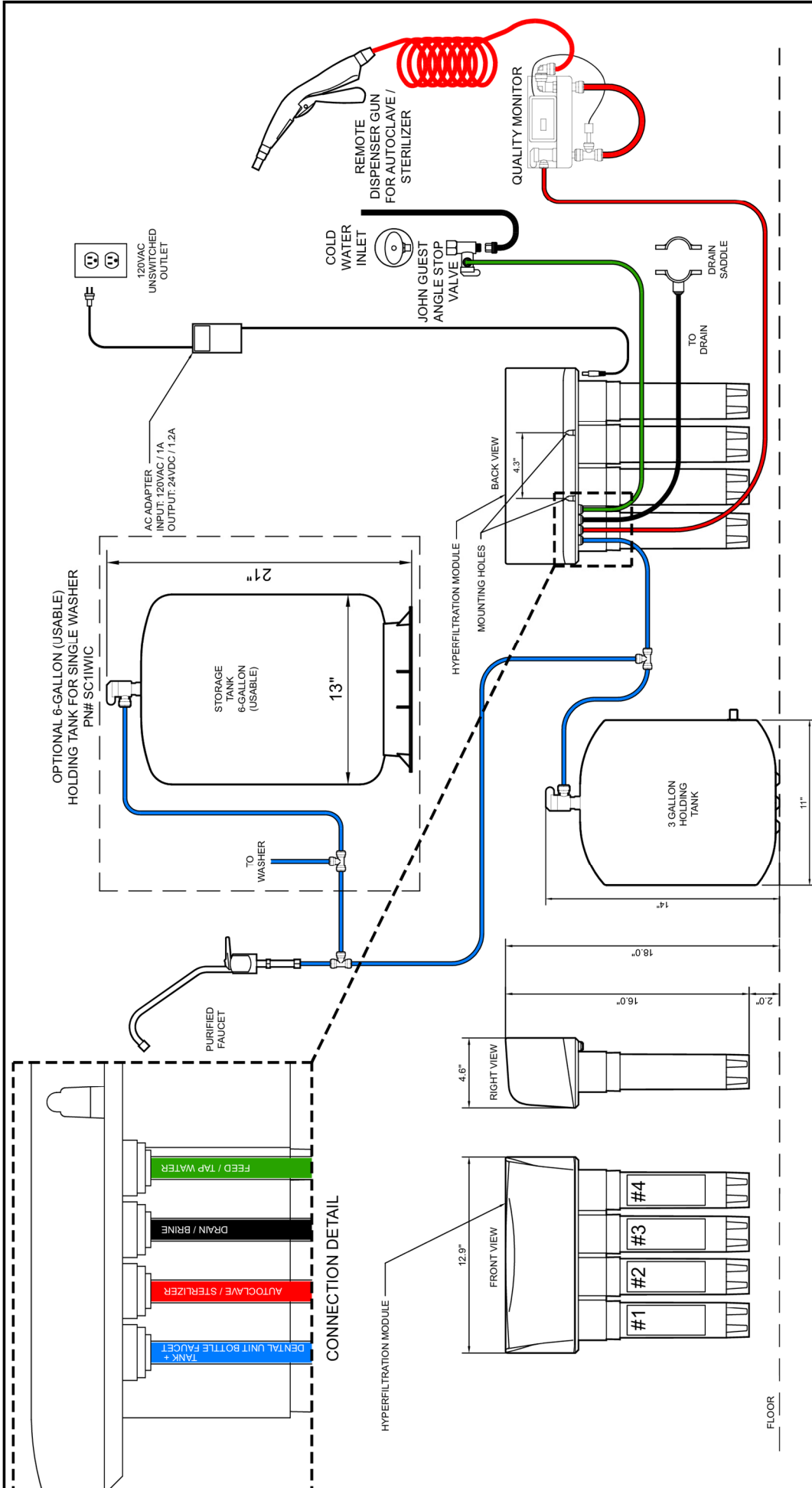

TITLE		SC-100		STERI-CENTER SYSTEM	
DATE	08/01/22	SIZE	B	SHEET	1 OF 1
SCALE	DO NOT SCALE	REV	-	PART NO.	

Note: Do not use with water that is microbiologically unsafe or of unknown quality. This device is only intended to be use as prescribed in the operation manual any miss-use, abuse improper installation or maintenance voids any and all warranties and claims.